



# Nucleoside Triphosphates (NTPs) Emprove® Expert


mRNA therapies represent a revolutionary approach to therapeutics such as vaccines, cancer treatments, and immunotherapies, offering novel ways to treat complex diseases. As a result, global interest in mRNA manufacturing has grown, while regulatory expectations have become increasingly stringent. In this environment, manufacturers require raw materials suitable for pharmaceutical use, ensuring the highest levels of quality and purity for successful mRNA production.

In response, we have developed the high-purity NTP Emprove® Expert portfolio, compliant with IPEC-PQG GMP standards, to ensure batch-to-batch consistency in your mRNA manufacturing processes. These critical raw materials are accompanied by comprehensive Emprove® dossiers, providing essential regulatory support to navigate the complexities of the biopharmaceutical industry while ensuring compliance and accelerating time to market.

The NTP Emprove® Expert portfolio includes the following products:

- Adenosine 5'-Triphosphate (ATP), 100 mM Sodium Solution EMPROVE® EXPERT
- Cytidine 5'-Triphosphate (CTP), 100 mM Sodium Solution EMPROVE® EXPERT
- Guanosine 5'-Triphosphate (GTP), 100 mM Sodium Solution EMPROVE® EXPERT
- Uridine 5'-Triphosphate (UTP), 100 mM Sodium Solution EMPROVE® EXPERT
- N1-Methyl-pseudouridine 5'-Triphosphate (N1-Me-pUTP), 100 mM Sodium Solution EMPROVE® EXPERT\*

The NTP Emprove® Expert products are provided as **100 mM sodium salt solutions**.

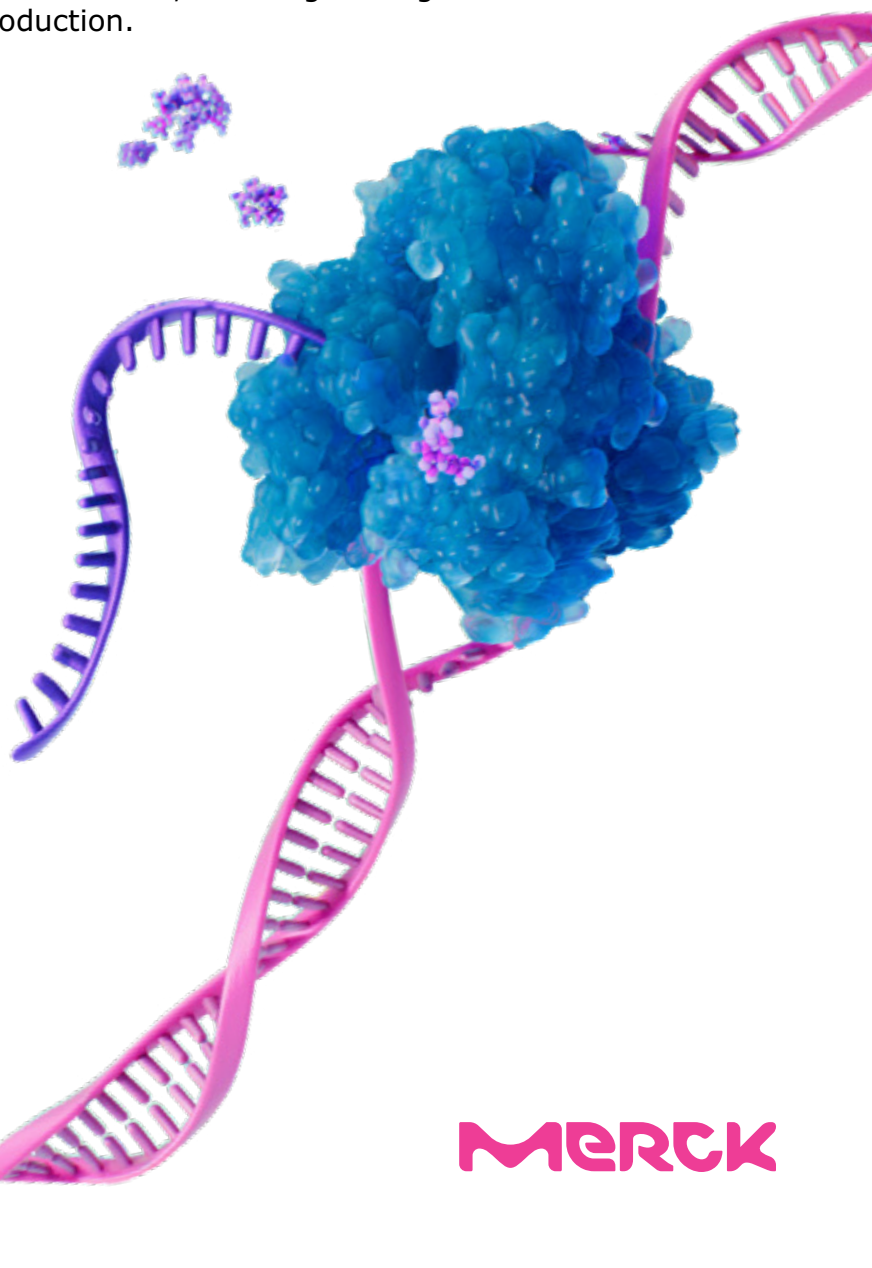

## Key Benefits

- **Improved batch-to-batch consistency** through IPEC-PQG GMP compliant manufacturing with validated QC methods, manufacturing processes, and stability studies
- **Animal origin-free (AOF)** to minimize contamination risks and support safety of high-risk applications.
- **High purity (≥99%) Emprove® Expert products** with specifications that meet the high demands of the biopharma industry:
  - **Specified low microbial and endotoxin levels**, suitable for high-risk applications.
  - **Specified absence of nuclease activity** to minimize degradation risk.
  - Extended impurity testing including specified low levels of **residual solvents**.
- **Backed by our Emprove® Program**, provides comprehensive regulatory, technical, and supply information to facilitate your risk assessment, material qualification, and process optimization efforts.

## Product Application

NTPs are crucial for producing mRNA-based therapeutics and vaccines. The four **natural/non-modified NTPs (ATP, CTP, GTP, UTP)** serve as essential building blocks for synthesizing mRNA through the in vitro

transcription (IVT) process, where T7 RNA polymerases incorporate them into the growing mRNA strand.

Non-modified NTPs can be exchanged by modified NTPs.<sup>[1,2]</sup>

## Product Specifications

| Name                    | ATP                                                                           | CTP                                                                           | GTP                                                                          | UTP                                                                          | N1-Me-pUTP                                                                    |
|-------------------------|-------------------------------------------------------------------------------|-------------------------------------------------------------------------------|------------------------------------------------------------------------------|------------------------------------------------------------------------------|-------------------------------------------------------------------------------|
| <b>Chemical Name</b>    | Adenosine Triphosphate                                                        | Cytidine Triphosphate                                                         | Guanosine Triphosphate                                                       | Uridine Triphosphate                                                         | N1-Methyl-Pseudouridine Triphosphate*                                         |
| <b>CAS Number</b>       | 56-65-5                                                                       | 86-01-1                                                                       | 65-47-4                                                                      | 63-39-8                                                                      | 1428903-59-6                                                                  |
| <b>Molecular Weight</b> | 507.18 g/mol                                                                  | 523.18 g/mol                                                                  | 483.18 g/mol                                                                 | 484.18 g/mol                                                                 | 498.17 g/mol                                                                  |
| <b>Chemical Formula</b> | C <sub>10</sub> H <sub>12</sub> N <sub>5</sub> O <sub>13</sub> P <sub>3</sub> | C <sub>10</sub> H <sub>15</sub> N <sub>5</sub> O <sub>14</sub> P <sub>3</sub> | C <sub>9</sub> H <sub>14</sub> N <sub>5</sub> O <sub>12</sub> P <sub>2</sub> | C <sub>9</sub> H <sub>13</sub> N <sub>2</sub> O <sub>12</sub> P <sub>2</sub> | C <sub>10</sub> H <sub>13</sub> N <sub>2</sub> O <sub>12</sub> P <sub>2</sub> |

*Molecular weight, CAS, Chemical formula is for free acid.*

## Physicochemical Information

- All NTPs are provided in a **sodium salt form**, dissolved in an **aqueous solution** at a concentration of **100 ± 5 mM and pH 7 ± 0.1**.
- High purity (≥99%)
- Storage: -20 °C
- Shelf-life: 1 year; extension to 2 years after completion of stability studies

## The Emprove® Program

### The Smart Way to Master Compliance and Control

Ensuring the compliance of your pharma and biopharma products involves the compilation of a vast amount of data, which can be time-consuming and resource-intensive. To help facilitate and accelerate your risk assessment continuum, we have developed the Emprove® Program, offering convenient access to reliable information for a broad portfolio of products.

Each product in the portfolio is complemented with different types of dossiers: Material Qualification Dossier, Quality Management Dossier, and Operational Excellence Dossier. They provide information on the

manufacturing process, stability data, elemental impurity information, product quality reports, analytical procedures, and much more.

The dossiers can be accessed online in our new Emprove® Suite, our information-as-a-service digital platform. A subscription can help you stay current. In addition to viewing and downloading dossiers, you can also receive notification updates on changes to documents, as well as generate metrics and reports.

For more information, please visit:  
[SigmaAldrich.com/emprove-program](https://SigmaAldrich.com/emprove-program)

## Ordering Information

| Product Description                                                                                 | Material Number              | Pack Size (mL) |
|-----------------------------------------------------------------------------------------------------|------------------------------|----------------|
| Adenosine 5'-Triphosphate (ATP),<br>100 mM Sodium Solution EMPROVE® EXPERT                          | 137201.0002                  | 1              |
|                                                                                                     | 137201.0010                  | 10             |
|                                                                                                     | 137201.0050                  | 50             |
|                                                                                                     | 137201.0200 (on demand only) | 200            |
| Guanosine 5'-Triphosphate (GTP),<br>100 mM Sodium Solution EMPROVE® EXPERT                          | 137202.0002                  | 1              |
|                                                                                                     | 137202.0010                  | 10             |
|                                                                                                     | 137202.0050                  | 50             |
|                                                                                                     | 137202.0200 (on demand only) | 200            |
| Cytidine 5'-Triphosphate (CTP),<br>100 mM Sodium Solution EMPROVE® EXPERT                           | 137203.0002                  | 1              |
|                                                                                                     | 137203.0010                  | 10             |
|                                                                                                     | 137203.0050                  | 50             |
|                                                                                                     | 137203.0200 (on demand only) | 200            |
| Uridine 5'-Triphosphate (UTP),<br>100 mM Sodium Solution EMPROVE® EXPERT                            | 137204.0002                  | 1              |
|                                                                                                     | 137204.0010                  | 10             |
|                                                                                                     | 137204.0050                  | 50             |
|                                                                                                     | 137204.0200 (on demand only) | 200            |
| N1-Methyl-pseudouridine 5'-Triphosphate<br>(N1-Me-pUTP), 100 mM Sodium Solution<br>EMPROVE® EXPERT* | 137206.0002                  | 1              |
|                                                                                                     | 137206.0010                  | 10             |
|                                                                                                     | 137206.0050                  | 50             |
|                                                                                                     | 137206.0200 (on demand only) | 200            |

## References

1. Karikó, K., et al. (2021). "The Critical Contribution of Pseudouridine to mRNA COVID-19 Vaccines." Nat Rev Immunol.
2. Li, B., et al. (2016). "Effects of Chemically Modified Messenger RNA on Protein Expression." Bioconjug Chem.

\* Customer shall at all times be solely responsible for obtaining all necessary licenses and rights to use this product for its intended application and any compliance with existing laws and regulations. Customer understands that the seller is not liable for any patent infringement claims related to Purchaser's use of this product.

We provide information and advice to our customers to the best of our knowledge and ability, but without obligation or liability.

Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

For additional information, please visit  
[SigmaAldrich.com /NTPs](https://www.sigmaaldrich.com/NTPs)

To place an order or receive technical assistance, please visit  
[SigmaAldrich.com/support/customer-support](https://www.sigmaaldrich.com/support/customer-support)

Merck KGaA  
Frankfurter Strasse 250  
64293 Darmstadt  
Germany

