

Data Sheet

Mouse OP9-DL4 Notch Ligand Cell Line (-GFP)

SCC495

Pack Size: $\geq 1 \times 10^6$ viable cells/vial

Store in liquid nitrogen.

FOR RESEARCH USE ONLY

Not for use in diagnostic procedures. Not for human or animal consumption.

Background

OP9-DL4 is a stable bone marrow derived stromal cell line expressing the Notch ligand- Delta like 4 (DL4) ectopically.

Notch signaling controls multiple cell-fate decisions. Four Notch receptors (Notch 1-Notch 4) have been identified in mammals. They interact with Jagged or Delta-like (DL) family members of Notch ligands. This is followed by cleavage of the intracellular domain of Notch and its subsequent translocation into the nucleus, where it binds with the transcription factors and activates transcription of various downstream target genes.¹

Notch signaling is essential in early T-cell lineage commitments. Bone marrow progenitor cells expressing constitutively active Notch have been shown to develop into CD4 and CD8 double positive T cells rather than B-cells. They also play a key role in development of CD4⁺ and CD8⁺ single positive cells from CD4 and CD8 double positive precursor thymocytes. Also, this signaling has been shown to promote development of T-cells with $\alpha\beta$ + TCRs. Thus, Notch signaling is critical for the earliest stages of T-cell commitment.²

Of the 4 Notch receptors, Notch1 (N1) receptor signaling has been shown to be sufficient for T cell development. The delta like ligands are physiologically relevant N1 ligands. DL1 interacts with both N1 and Notch 2 (N2) to induce T-cell lineage commitment. However, DL4 interacts specifically only with N1 and supports T cell commitment and maturation both *in vitro* and *in vivo*. Also, results from previous binding studies show that binding between DL1 and N1 is weak as opposed to the stronger binding between DL4 and N1.³ Of the DL1 and DL4 delta like ligands, DL4 exhibits greater capacity to activate Notch pathway in hematopoietic progenitor cells.⁴

The OP9 bone marrow stromal cells that support the differentiation of hematopoietic progenitor cells (HPCs) into multiple lineages including B-cells, failed to support differentiation of HPCs to T-cells. This is mainly because they do not express Delta-like ligands. Therefore, OP9 cells were transduced with pMiGR retroviral vector engineered to express the Delta-like-4. These OP9-DL4 cells support the differentiation of embryonic or hematopoietic stem cells from fetal liver or bone marrow into T lymphocytes.

The OP9-DL4 Cell line (SCC495) has had the -GFP coding region removed from the original pMiGR plasmid, resulting in no -GFP expression for this cell line.

Source

Parental OP9-DL4 cells were genetically modified from OP9 stromal cells derived from mouse bone marrow.

Short Tandem Repeat

M18-3:16	M1-2: 17	M8-1: 17	M11-2: 16	MX-1: 28
M4-2:19.3, 20.3	M7-1: 26.2	M2-1: 15, 16	M17-2: 15	M13-1: 17
M6-7: 17	M1-1: 16, 17	M15-3: 22.3, 26.3	M12-1: 17	
M19-2: 13	M3-2: 13, 14	M6-4: 18	M5-5: 15	

Quality Control Testing

- The OP9-DL4 without -GFP murine bone marrow stromal cells are verified to be of mouse origin and negative for human, rat, Chinese hamster, Golden Syrian hamster, and nonhuman primate interspecies contamination, as assessed by a Contamination Clear panel by Charles River Animal Diagnostic Services.
- Cells tested negative for infectious diseases against a Mouse Essential CLEAR panel by Charles River Animal Diagnostic Services.
- Cells tested negative for mycoplasma.

Storage and Handling

The OP9-DL4 without -GFP cells should be stored in liquid nitrogen until use. The cells can be cultured for at least 10 passages after initial thawing without significantly affecting the cell marker expression and functionality.

Representative Data

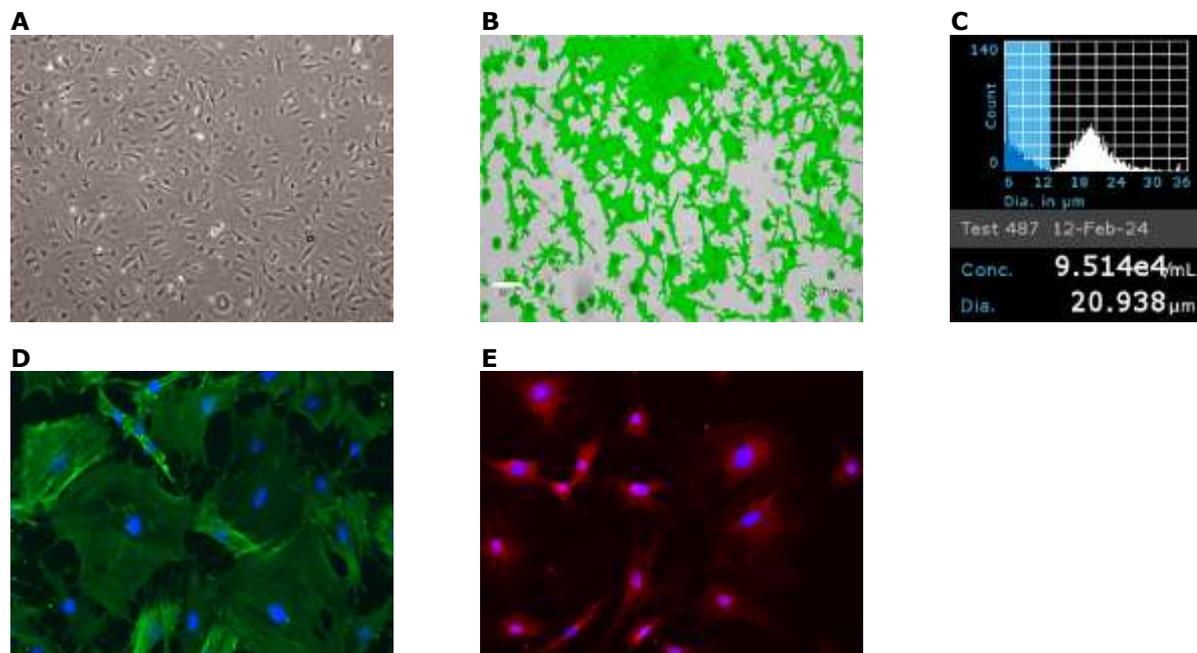


Figure 1. (A) Bright-field image of OP9-DL4 without -GFP cells one day after thaw in a T225 flask (4X magnification). (B) Cell confluency was assessed throughout the culture using the Millicell® Digital Cell Imager (MDCI10000). (C) Cell counting was performed using the Scepter™ 3.0 Handheld Automated Cell Counter using 60 μm sensors (PHCC360KIT). (D) The OP9-DL4 without -GFP cells stained with Phalloidin-Atto-647N (65906). (E) OP9-DL4 without -GFP cells express DL4 protein (HPA023392).

Note: Product catalog numbers indicated in () can be purchased at [SigmaAldrich.com](https://www.sigmaaldrich.com) unless otherwise stated.

Protocols

Thawing the Cells

1. Do not thaw the cells until the recommended medium is on hand. Cells can grow on standard tissue cultureware surfaces without any additional coating.
2. The OP9-DL4 without -GFP cells are thawed and expanded in OP9-DL4 Expansion Medium comprising of MEM-Eagle, alpha modification (M4526) containing 5% FBS (ES-009-B), 2 mM L-Glutamine (G7513) and Penicillin/Streptomycin (P4333) (optional).
3. Remove the vial of frozen OP9-DL4 without -GFP cells from liquid nitrogen and incubate in a 37 °C water bath. Closely monitor until the cells are completely thawed. Maximum cell viability is dependent on the rapid and complete thawing of frozen cells.
IMPORTANT: Do not vortex the cells.
4. As soon as the cells are completely thawed, disinfect the outside of the vial with 70% ethanol. Proceed immediately to the next step.
5. In a laminar flow hood, use a 1-2 mL pipette to transfer the cells to a sterile 15 mL conical tube. Be careful not to introduce any bubbles during the transfer process.
6. Using a 10 mL pipette, slowly add dropwise 9 mL of OP9-DL4 Expansion Medium (Step 1 above) to the 15 mL conical tube.
IMPORTANT: Do not add the entire volume of media all at once to the cells. This may result in decreased cell viability due to osmotic shock.
7. Gently mix the cell suspension by slowly pipetting up and down twice. Be careful not to introduce any bubbles.
IMPORTANT: Do not vortex the cells.
8. Centrifuge the tube at 300 x *g* for 2-3 minutes to pellet the cells.
9. Decant as much of the supernatant as possible. Steps 5-8 are necessary to remove residual cryopreservative (DMSO).
10. Resuspend the cells in 35 mL of OP9-DL4 Expansion Medium.
11. Transfer the cell mixture to a T225 tissue culture flask.
12. Incubate the cells at 37 °C in a humidified incubator with 5% CO₂.

Subculturing the Cells

1. The OP9-DL4 without -GFP cells can be passaged at ~80-85% confluency.
2. Carefully remove the medium from the tissue culture flask containing the 80-85% confluent layer of OP9-DL4 without -GFP cells.
3. Rinse the flask with 10 mL 1X sterile PBS (TMS-012-A). Aspirate after the rinse.
4. Apply 5-7 mL of pre-warmed Accutase® (A6964) and incubate in a 37 °C incubator for 5 minutes.
5. Inspect the flask and ensure the complete detachment of cells by gently tapping the side of the flask with the palm of your hand.
6. Add 5-7 mL of OP9-DL4 Expansion Medium to the plate.
7. Gently rotate the flask to mix the cell suspension. Transfer the dissociated cells to a 15 mL conical tube.
8. Centrifuge the tube at 300 x *g* for 3-5 minutes to pellet the cells.
9. Discard the supernatant, then loosen the cell pellet by tapping the tip of the tube with a finger.
10. Apply 2-5 mL of OP9-DL4 cell medium to the conical tube and resuspend the cells thoroughly. Large cell clumps may be broken up by gentle trituration.
IMPORTANT: Do not vortex the cells.
11. Count the number of cells using a hemocytometer or a Scepter™ 3.0 Handheld Automated Cell Counter.
12. Plate the cells to the desired density. Typical split ratio is 1:6.

Cryopreservation of the Cells

The OP9-DL4 without -GFP cells may be frozen in OP9-DL4 Expansion Medium supplemented with 10% DMSO using a Nalgene® slow freeze Mr. Frosty® container.

References

1. Schmitt TM and Zúñiga-Pflücker JC. 2002. Induction of T cell development from hematopoietic progenitor cells by delta-like-1 in vitro. *Immunity*. 17(6): 749–756.
2. de Pooter R and Zúñiga-Pflücker JC. 2007. T-cell potential and development in vitro: the OP9-DL1 approach. *Curr Opin Immunol*. 19(2):163–168.
3. Besseyrias V, Fiorini E, Strobl LJ, Zimmer-Strobl U, Dumortier A, Koch U, Arcangeli M-L, Ezine S, Macdonald HR, Radtke F. 2007. Hierarchy of Notch–Delta interactions promoting T cell lineage commitment and maturation. *J Exp Med*. 204(2):331–343.
4. Mohtashami M, Shah DK, Nakase H, Kianizad K, Petrie HT, Zúñiga-Pflücker JC. 2010. Direct comparison of Dll1- and Dll4-mediated Notch activation levels shows differential lymphomyeloid lineage commitment outcomes. *J Immunol*. 185(2):867–876.

Academic Use Agreement

Subject to local law

THIS PRODUCT MAY ONLY BE USED BY INDIVIDUALS EMPLOYED BY AN ACADEMIC INSTITUTION AND IS INTENDED SOLELY TO BE USED FOR ACADEMIC RESEARCH, WHICH IS FURTHER DEFINED BELOW. BY OPENING THIS PRODUCT, YOU ("PURCHASER") HEREBY REPRESENT THAT YOU HAVE THE RIGHT AND AUTHORITY TO LEGALLY BIND YOURSELF AND/OR YOUR EMPLOYER INSTITUTION, AS APPLICABLE, AND CONSENT TO BE LEGALLY BOUND BY THE TERMS OF THIS ACADEMIC USE AGREEMENT. IF YOU DO NOT AGREE TO COMPLY WITH THESE TERMS, YOU MAY NOT OPEN OR USE THE PRODUCT AND YOU MUST CALL CUSTOMER SERVICE (1-800-645-5476) TO ARRANGE TO RETURN THE PRODUCT FOR A REFUND.

"Product" means Mouse OP9-DL4 Notch Ligand Cell Line (-GFP) (SCC495).

"Academic Research" means any internal in vitro research use by individuals employed by an academic institution. Academic Research specifically excludes the following uses of whatever kind or nature:

- Re-engineering or copying the Product.
- Making derivatives, modifications, or functional equivalents of the Product
- Obtaining patents or other intellectual property rights claiming use of the Product
- Using the Product in the development, testing, or manufacture of a Commercial Product
- Using the Product as a component of a Commercial Product
- Reselling or licensing the Product.
- Using the Product in clinical or therapeutic applications including producing materials for clinical trials
- Administering the Product to humans
- Using the Product in collaboration with a commercial or non-academic entity

"Commercial Product" means any product intended for: (i) current or future sale; (ii) use in a fee-for-service; or (iii) any diagnostic, clinical, or therapeutic use.

Access to the Product is limited solely to those officers, employees, and students of PURCHASER's academic institution who need access to the Product to perform Academic Research. PURCHASER shall comply with all applicable laws in its use and handling of the Product and shall keep it under reasonably safe and secure conditions to prevent unauthorized use or access.

These use restrictions will remain in effect for as long as PURCHASER possesses the Product.

COMMERCIAL OR NON-ACADEMIC ENTITIES INTERESTED IN PURCHASING OR USING THE PRODUCT MUST CONTACT licensing@milliporesigma.com AND AGREE TO SEPARATE TERMS OF USE PRIOR TO USE OR PURCHASE.

Genetically Modified Organisms (GMO)

This product contains genetically modified organisms.
Este producto contiene organismos genéticamente modificados.
Questo prodotto contiene degli organismi geneticamente modificati.
Dieses Produkt enthält genetisch modifizierte Organismen.
Ce produit contient des organismes génétiquement modifiés.
Dit product bevat genetisch gewijzigde organismen.
Tämä tuote sisältää geneettisesti muutettuja organismeja.
Denna produkt innehåller genetiskt ändrade organismer.

Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The information in this document is subject to change without notice and should not be construed as a commitment by the manufacturing or selling entity, or an affiliate. We assume no responsibility for any errors that may appear in this document.

Technical Assistance

Visit the tech service page at SigmaAldrich.com/techservice.

Terms and Conditions of Sale

Warranty, use restrictions, and other conditions of sale may be found at SigmaAldrich.com/terms.

Contact Information

For the location of the office nearest you, go to SigmaAldrich.com/offices.

The life science business of Merck operates
as MilliporeSigma in the U.S. and Canada.

Merck Millicell, Scepter and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

© 2024 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

Document Template 20306518 Ver 6.0

23060975 Ver 1.0, Rev 02MAY2024, AV

The Merck logo is displayed in a bold, red, sans-serif font.