

Data Sheet

## Mouse EERL Lung Metastasis (MLM) Cell Line

### SCC625

Pack Size:  $\geq 1 \times 10^6$  viable cells/vial

Store in liquid nitrogen.

#### FOR RESEARCH USE ONLY

**Not for use in diagnostic procedures. Not for human or animal consumption.**

### Background

The head and neck cancers represent the sixth most common class of cancer.<sup>1,2</sup> Among them, incidents of the oropharyngeal squamous cell carcinoma (OPSCC) are on the rise, and a majority of them are caused by the human papilloma virus (HPV).<sup>3</sup> While successful management methods exist, the local recurrence and the nodal metastasis do result in mortality in a minority of the cases. Historically, insights into the HPV-induced transformation have been obtained in the context of the cervical cancer and the anogenital warts models, and, while much of it is relevant to OPSCC, it will be beneficial to have the systems in which the effects of HPV infection can be studied in the context specific to OPSCC.

The mEER cell line (SCC626), was established by the retroviral transduction of the HPV E6 gene, the HPV E7 gene, and the human hRas G12V gene into primary mouse oropharyngeal epithelial cells.<sup>4</sup> The resulting cell line was shown to be capable of forming poorly differentiated, aggressively metastatic tumors following the orthotopic injections in immune-competent mice. SCC627, (mEERL), the luciferase-expressing derivative of SCC626, was successfully used to optically track the metastasis through the lymphatic space to the lungs, which is anatomically relevant to the human oral cancers.<sup>4,5</sup>

The mEERL Lung Metastasis cell line (MLM, SCC625), is one of the five clones isolated from such lung metastasis of the SCC627 inoculation that formed after the apparent disease clearance following a treatment by cisplatin and radiation.<sup>6</sup> The mEERL Lung Metastasis cell line (SCC625), was shown to be inherently more resistant to the same treatments than its parent and showed ready tropism to the lungs upon re-injection into immune-competent mice, a higher rate of metastasis, and a higher number of metastatic tumors than its parent or the four other siblings isolated in parallel. While this clone has lost the luciferase expression of the parental SCC627, it retained the expression of the many of the so-called hallmark genes of OPSCC, including cytokeratin, E-cadherin, BRCA2, p16, and EGFR, and its spread has been tracked *in vivo* by such markers.<sup>6</sup>

### Source

The mEERL (SCC625) cell line is a clone isolated from the lung metastasis of the mEERL cell line (SCC627) in a treatment-failed male C57BL/6 mouse.<sup>5,6</sup> SCC627 is a luciferase-expressing derivative of the mEER cell line (SCC626). The mEERL cell line (SCC626) was established by the retroviral transduction of the HPV16 E6 gene, HPV16 E7 gene, and the human hRas gene into the primary tonsillar epithelial cells isolated from a male C57BL/6 mouse.<sup>4</sup>

## Short Tandem Repeat

M18-3:	16,17	M1-2:	18	M8-1:	16	M11-2:	16,17	MX-1:	28
M4-2:	20.3	M7-1:	26.2	M2-1:	16	M17-2:	15	M13-1:	17
M6-7:	15	M1-1:	15,16	M15-3:	22.3	M12-1:	17		
M19-2:	13	M3-2:	15	M6-4:	18	M5-5:	17		


## Quality Control Testing

- The MLM cell line is verified to be of mouse origin and negative for human, rat, Chinese hamster, Golden Syrian hamster, and non-human primate interspecies contamination, as assessed by a Contamination Clear panel by Charles River Animal Diagnostic Services.
- Cells tested negative for infectious diseases against a Mouse Essential CLEAR panel by Charles River Animal Diagnostic Services.
- Cells tested negative for mycoplasma.


## Storage and Handling

The MLM cells should be stored in liquid nitrogen until use. The cells can be cultured for at least 10 passages after initial thawing without significantly affecting the cell marker expression and functionality.


## Representative Data

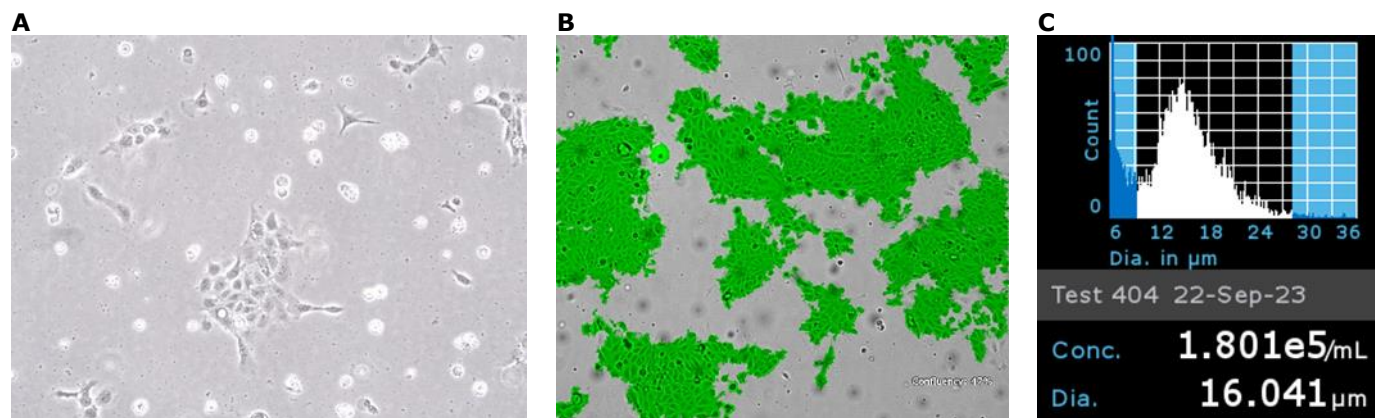

**Figure 1.** (A) Bright-field images of MLM cells a day after thaw in a T175 flask. (B) Cell confluency was assessed throughout the culture using the Millicell® Digital Cell Imager (MDCI10000). (C) Cell counting was performed using Scepter™ 3.0 Handheld Automated Cell Counter using 60 μm sensor tips (PHCC360KIT).

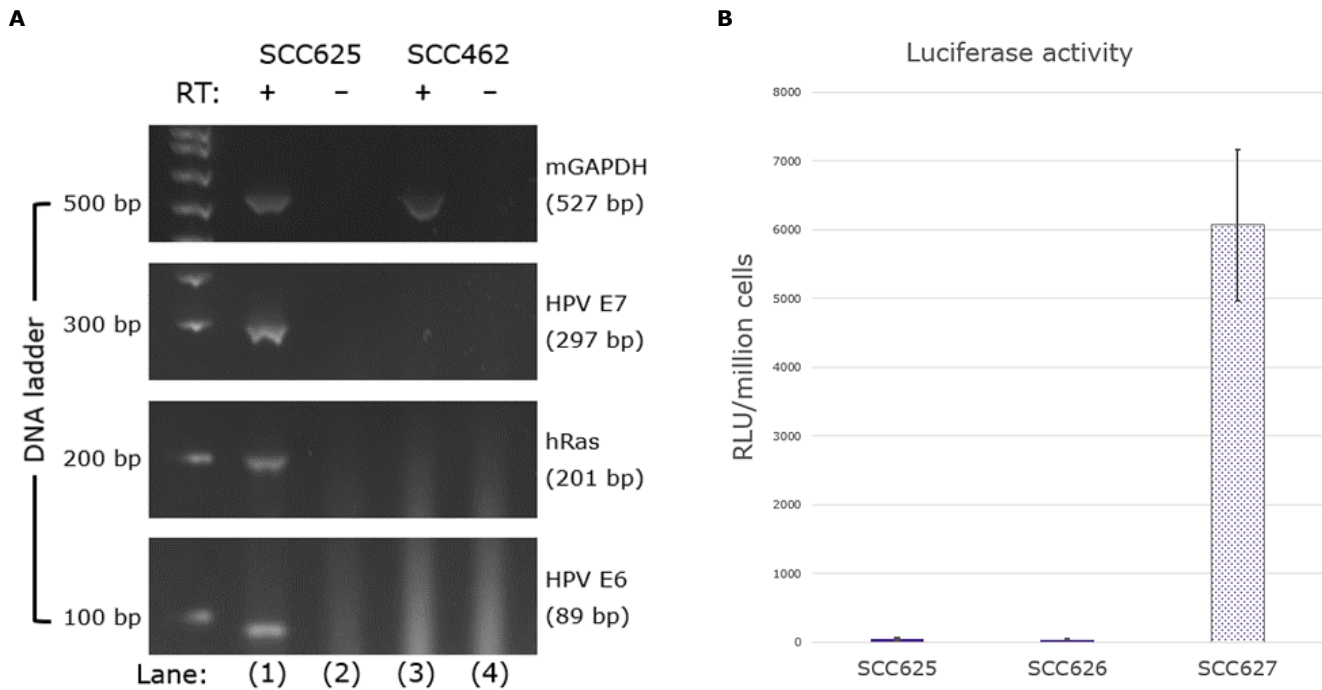

**Figure 2.** Transgene expression in SCC625, the mouse EERL Lung Metastasis (MLM) cell line. **(A)** End-point RT-PCR with primers specific for HPV16 E6 (expected product size 89 bp), HPV16 E7 (expected 297 bp), and human hRas (expected 201 bp) successfully detected these transcripts, when the total RNA prepared from SCC625 cells was pre-incubated with reverse transcriptase (Lane 1) but not without the RT enzyme (Lane 2). These transgenes were not detected in an untransfected (and unrelated) mouse cell line, with or without RT (SCC462, Ret-melanoma Sorted cell line, Lanes 3 and 4). Mouse GAPDH (expected 527 bp) was used as control. The TaqMan RNA-to-CT 1-Step Kit (Applied Biosystems 4392938) was used for the reactions. **(B)** Expression of the luciferase gene in the parental SCC627, the mEERL cell line, was reportedly lost in the lung metastasis clone from which SCC625 was derived, and that result is also confirmed here.<sup>6</sup> Luciferase activity was measured using the Promega Luciferase Assay System (Fisher Scientific, PR-E 1500). Each bar represents N=3 +/- standard deviation.

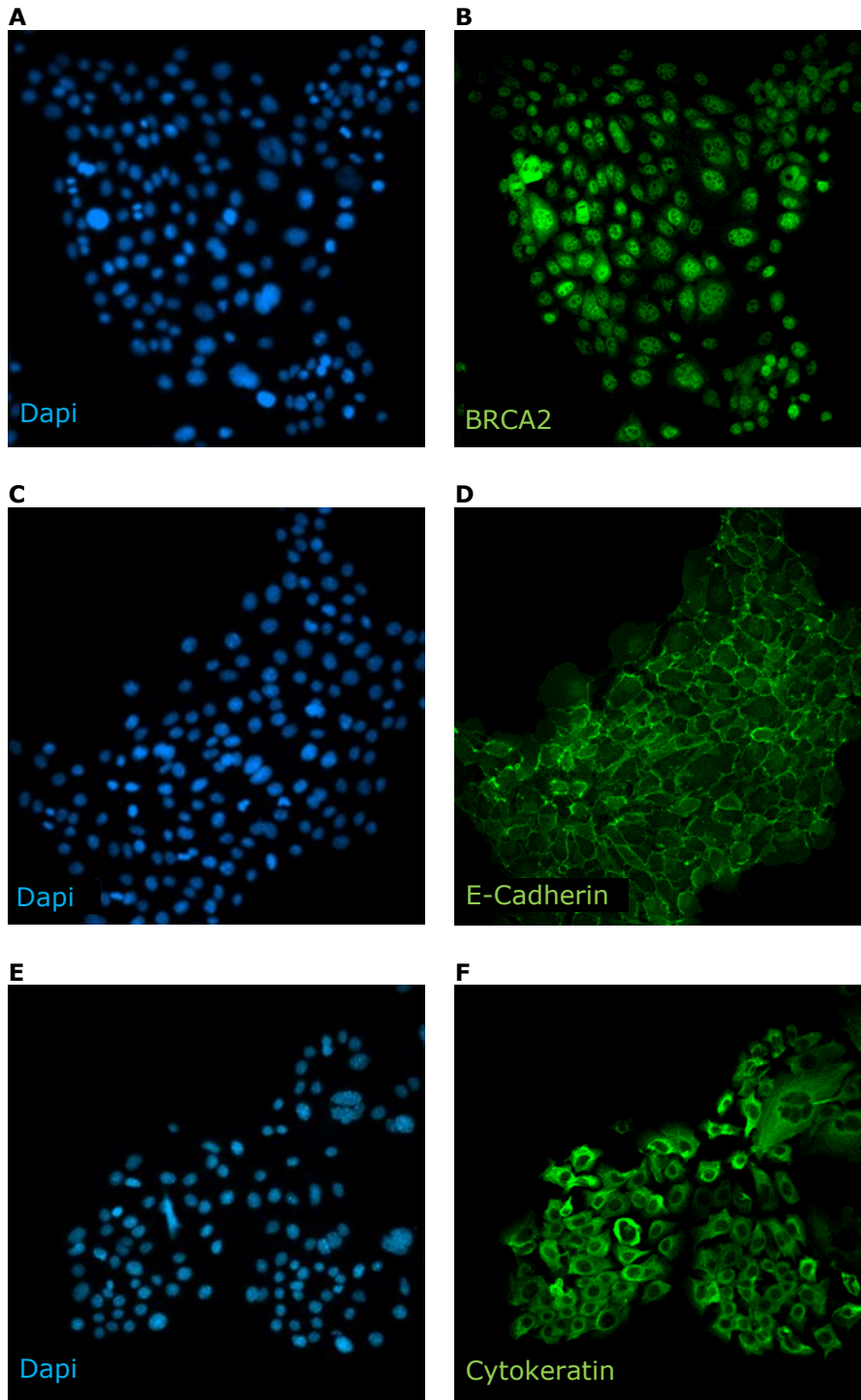

**Figure 3.** (A, C, E) MLM cells express hallmark genes of oropharyngeal squamous cell carcinoma. MLM cells labeled with DAPI (MBD0015), (B) anti-BRCA2 antibodies (HPA026815), (D) anti-E-Cadherin antibodies (SAB5700789), and (F) anti-cytokeratin antibodies (C2931).

## Protocols

The MLM cells proliferate rapidly. We recommend thawing in a T175 flask and use of trypsin to detach the cells.

### Thawing the Cells

1. Do not thaw the cells until the recommended medium is on hand. Cells can grow on standard tissue cultureware surfaces without any additional coating.  
The MLM cells are thawed and expanded in MLM Expansion Medium comprising of DMEM/F12 (DF-042-B) containing 10% FBS (ES-009-B), 0.5 µg/mL hydrocortisone (H0888), 5 µg/mL human apo-Transferrin (T4382), 5 µg/mL human insulin (I9278), 5 ng/mL human Epidermal Growth Factor (SRP3027), 1.36 ng/mL triiodothyronine (T6397), with optional Penicillin/Streptomycin (P4333).
2. Remove the vial of frozen cells from liquid nitrogen and incubate in a 37 °C water bath. Closely monitor until the cells are completely thawed. Maximum cell viability is dependent on the rapid and complete thawing of frozen cells.  
**Important:** Do not vortex the cells.
3. As soon as the cells are completely thawed, disinfect the outside of the vial with 70% ethanol. Proceed immediately to the next step.
4. In a laminar flow hood, use a 1- or 2- mL pipette to transfer the cells to a sterile 15 mL conical tube. Be careful not to introduce any bubbles during the transfer process.
5. Using a 10 mL pipette, slowly add dropwise 4 mL of MLM Expansion Medium (medium composition in Step 1) to the 15 mL conical tube.  
**Important:** Do not add the entire volume of media all at once to the cells. This may result in decreased cell viability due to osmotic shock.
6. Gently mix the cell suspension by slowly pipetting up and down twice. Be careful not to introduce any bubbles.  
**Important:** Do not vortex the cells.
7. Centrifuge the tube at 300 x *g* for 5 minutes to pellet the cells.
8. Decant as much of the supernatant as possible. Steps 5-8 are necessary to remove residual cryopreservative (DMSO).
9. Resuspend the cells in 40 mL of MLM Expansion Medium.
10. Transfer the cell mixture to a T175 tissue culture flask.
11. Incubate the cells at 37 °C in a humidified incubator with 5% CO<sub>2</sub>.

### Subculturing the Cells

1. The MLM cells can be passaged at ~80-85% confluency.
2. Carefully remove the medium from the T175 tissue culture flask containing the 80-85% confluent layer of cells.
3. Rinse the flask with 30-40 mL 1X sterile PBS (TMS-012-A). Aspirate after the rinse. Repeat this wash step one more time.
4. Apply 10 mL of 0.25% Trypsin-EDTA pre warmed to 37 °C (SM-2003-C) and incubate in a 37 °C incubator for 5-7 minutes.
5. Inspect the flask and ensure the complete detachment of cells by gently tapping the side of the flask with the palm of your hand.
6. Add 10-15 mL of MLM Expansion Medium to the flask.
7. Gently rotate the flask to mix the cell suspension. Transfer the dissociated cells to a 50 mL conical tube.
8. Centrifuge the tube at 300 x *g* for 5 minutes to pellet the cells.
9. Discard the supernatant, then loosen the cell pellet by tapping the tip of the tube with a finger.
10. Apply 2-5 mL of MLM expansion medium to the conical tube and resuspend the cells thoroughly. Large cell clumps may be broken up by gentle trituration.  
**Important:** Do not vortex the cells.
11. Count the number of cells using a hemocytometer or a Scepter™ 3.0 Handheld Automated Cell Counter.
12. Plate the cells to the desired density. Typical split ratio is 1:10.

## Cryopreservation of the Cells

The MLM cells may be frozen in 90% fetal bovine serum (FBS) with 10% DMSO using a Nalgene® slow freeze Mr. Frosty® container.

## References

1. Kademani D. 2007. Oral cancer. Mayo Clin Proc. 82(7): 878-87.
2. Rogers SN, Brown JS, Woolgar JA, Lowe D, Magennis P, Shaw RJ, Sutton D, Errington D, Vaughan D. 2009. Survival following primary surgery for oral cancer. Oral Oncol. 45(3):201-11.
3. Pytynia KB, Dahlstrom KR, Sturgis EM. 2014. Epidemiology of HPV-associated oropharyngeal cancer. Oral Oncol. 50(5): 380-6.
4. Hoover AC, Spanos WC, Harris GF, Anderson ME, Klingelutz AJ, Lee JH. 2007. The role of human papillomavirus 16 E6 in anchorage-independent and invasive growth of mouse tonsil epithelium. Arch Otolaryngol Head Neck Surg. 133(5): 495-502.
5. Williams R, Lee DW, Elzey BD, Anderson ME, Hostager BS, Lee JH. 2009. Preclinical models of HPV+ and HPV- HNSCC in mice: an immune clearance of HPV+ HNSCC. Head Neck. 31(7): 911-8.
6. Vermeer DW, Coppock JD, Zeng E, Lee KM, Spanos WC, Onken MD, Uppaluri R, Lee JH, Vermeer PD. 2016. Metastatic model of HPV+ oropharyngeal squamous cell carcinoma demonstrates heterogeneity in tumor metastasis. Oncotarget. 7(17): 24194-207.

## Academic Use Agreement

Subject to local law

THIS PRODUCT MAY ONLY BE USED BY INDIVIDUALS EMPLOYED BY AN ACADEMIC INSTITUTION AND IS INTENDED SOLELY TO BE USED FOR ACADEMIC RESEARCH, WHICH IS FURTHER DEFINED BELOW. BY OPENING THIS PRODUCT, YOU ("PURCHASER") HEREBY REPRESENT THAT YOU HAVE THE RIGHT AND AUTHORITY TO LEGALLY BIND YOURSELF AND/OR YOUR EMPLOYER INSTITUTION, AS APPLICABLE, AND CONSENT TO BE LEGALLY BOUND BY THE TERMS OF THIS ACADEMIC USE AGREEMENT. IF YOU DO NOT AGREE TO COMPLY WITH THESE TERMS, YOU MAY NOT OPEN OR USE THE PRODUCT AND YOU MUST CALL CUSTOMER SERVICE (1-800-645-5476) TO ARRANGE TO RETURN THE PRODUCT FOR A REFUND.

"Product" means Mouse EERL Lung Metastasis (MLM) Cell Line (SCC625).

"Academic Research" means any internal in vitro research use by individuals employed by an academic institution. Academic Research specifically excludes the following uses of whatever kind or nature:

- Re-engineering or copying the Product
- Making derivatives, modifications, or functional equivalents of the Product
- Obtaining patents or other intellectual property rights claiming use of the Product
- Using the Product in the development, testing, or manufacture of a Commercial Product
- Using the Product as a component of a Commercial Product
- Reselling or licensing the Product
- Using the Product in clinical or therapeutic applications including producing materials for clinical trials
- Administering the Product to humans
- Using the Product in collaboration with a commercial or non-academic entity

"Commercial Product" means any product intended for: (i) current or future sale; (ii) use in a fee-for-service; or (iii) any diagnostic, clinical, or therapeutic use.

Access to the Product is limited solely to those officers, employees, and students of PURCHASER's academic institution who need access to the Product to perform Academic Research. PURCHASER shall comply with all applicable laws in its use and handling of the Product and shall keep it under reasonably safe and secure conditions to prevent unauthorized use or access.

These use restrictions will remain in effect for as long as PURCHASER possesses the Product.

COMMERCIAL OR NON-ACADEMIC ENTITIES INTERESTED IN PURCHASING OR USING THE PRODUCT MUST CONTACT [licensing@milliporesigma.com](mailto:licensing@milliporesigma.com) AND AGREE TO SEPARATE TERMS OF USE PRIOR TO USE OR PURCHASE.

## Genetically Modified Organisms (GMO)

This product contains genetically modified organisms.  
Este producto contiene organismos genéticamente modificados.  
Questo prodotto contiene degli organismi geneticamente modificati.  
Dieses Produkt enthält genetisch modifizierte Organismen.  
Ce produit contient des organismes génétiquement modifiés.  
Dit product bevat genetisch gewijzigde organismen.  
Tämä tuote sisältää geneettisesti muutettuja organismeja.  
Denna produkt innehåller genetiskt ändrade organismer.



---

## Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The information in this document is subject to change without notice and should not be construed as a commitment by the manufacturing or selling entity, or an affiliate. We assume no responsibility for any errors that may appear in this document.

## Technical Assistance

Visit the tech service page at [SigmaAldrich.com/techservice](https://SigmaAldrich.com/techservice).

## Terms and Conditions of Sale

Warranty, use restrictions, and other conditions of sale may be found at [SigmaAldrich.com/terms](https://SigmaAldrich.com/terms).

## Contact Information

For the location of the office nearest you, go to [SigmaAldrich.com/offices](https://SigmaAldrich.com/offices).

The life science business of Merck operates  
as MilliporeSigma in the U.S. and Canada.

Merck Millicell, Scepter and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

© 2024 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

Document Template 20306518 Ver 6.0

23060982 Ver 1.0, Rev 06MAY2024, DP

