# Millipore<sub>®</sub>

Preparation, Separation, Filtration & Monitoring Products





## FLASH® RT Allergen Indicator Protein Test


Cleaning Verification for Protein Removal and Allergen Control Programs.

FLASH® RT is a total protein visual test that rapidly detects protein residues left on food contact surfaces after cleaning. Protein is a difficult food residue to remove. Most allergens are proteins, so quick verification of surface hygiene helps minimize the risk of cross-contact to allergen-free products. Using FLASH® RT regularly also helps reduce the opportunity for biological contamination.

As part of a HACCP allergen management program, FLASH® RT supports process verification requirements that ensure cleaning methods, independently validated to effectively remove allergens, are consistently applied.

Used together with the MVP ICON® for ATP testing, both can be incorporated into validated cleaning procedures (SSOPs) as part of a more comprehensive HACCP management program.

To find out more, visit: SigmaAldrich.com

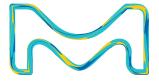

### A Perfect Complement to ATP Testing


#### **Colorimetric Assay. No Equipment Needed.**

FLASH® RT is a self-contained sampling device containing a pre-moistened tip and all necessary reagents. Testing is done at ambient room temperature incubation and requires no equipment.

#### Manage Risks. Be Prepared. Everytime.

- Detects total protein including "Big 8" food allergens
- Rapid pass / fail result within 10 minutes
- Detects down to 20 µg total protein
- Simple online training available 24/7 with certification for audit compliance

#### 3 Step Process. Simple.

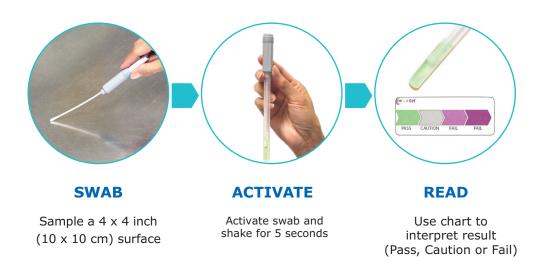

To place an order or receive technical assistance
In the U.S. and Canada, call toll-free 1-800-245-0113
For other countries across Europe and the world, please visit: SigmaAldrich.com/offices
For Technical Service, please visit: SigmaAldrich.com/techservice

Merck KGaA Frankfurter Strasse 250 64293 Darmstadt, Germany

SigmaAldrich.com

