

30574 Atto 488 Biotin


Application

Atto 488 is a new hydrophilic fluorescent label with excellent water solubility. The dye exhibits strong absorption, high fluorescence quantum yield and exceptional thermal and photo-stability. Thus Atto 488 is highly suitable for single-molecule detection applications and high-resolution microscopy such as PALM, dSTORM, STED etc. Additionally the dye highly qualifies to be applied in flow cytometry (FACS), fluorescence in-situ hybridization (FISH) and many more. The fluorescence is excited most efficiently in the range 480 - 515 nm. A suitable source of excitation is the 488 nm line of the Argon-Ion laser.


Biotin conjugates can be used in applications like ELISA or immuno-histochemistry, in situ hybridization, flow cytometry and others, to identify streptavidin, avidin, or extravidin-conjugates.

Product Description

MW	1191 g/mol
λ_{abs}	500 nm
ϵ_{max}	$9.0 \times 10^4 \text{ M}^{-1} \text{ cm}^{-1}$
λ_{fl}	520 nm
η_{fl}	80 %
τ_{fl}	4.1 ns
CF ₂₆₀	0.22
CF ₂₈₀	0.09

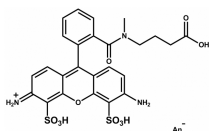

Optical data of the carboxy derivative (in aqueous solution)

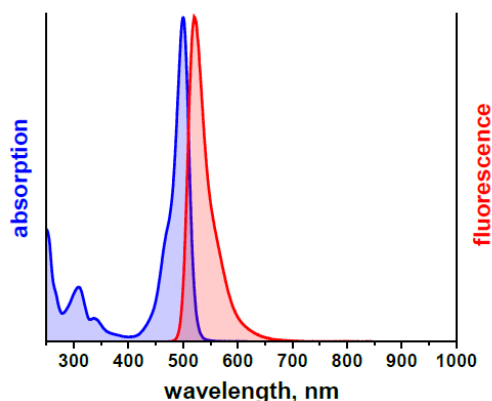

Storage: Store at -20°C and protected from light.

Precautions and Disclaimer

This product is for R&D use only, not for drug, household, or other uses. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.

The vibrant M and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates.

Detailed information on trademarks is available via publicly accessible resources.

© 2018 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

The life science business of Merck KGaA, Darmstadt, Germany operates as MilliporeSigma in the US and Canada.

