


93543**N-Hydroxysuccinimide (NHS) activated MegaStokes dye 605****Molecular formula :**C₂₄H₂₈F₆N₃O₄P**Molecular weight (g/mol) :**

567.46

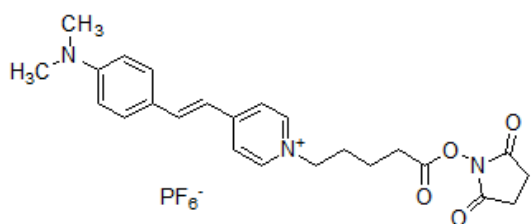

The MegaStokes dyes for „click“ biological labeling, licensed by LuminoChem are virtually covering the whole visible spectrum reaching the infra-red regime. Besides dyes that are being capable of participating in classical copper catalyzed 1,3-dipolar cycloaddition reaction with the counterparting function, we also provide dyes containing cyclooctyne moiety an alkyne derivative that enables copper free clicking to azides.

These dyes are noteworthy for their large Stokes-shift (120 or 125 nm). Dyes for this kind are exceptionally useful in fluorescence resonance electron transfer (FRET) applications as they do not interfere with the spectral bands of the second fluorophore a common problem in FRET technology. Azido sugars have been used for click-labeling of surface glycoproteins, recently. We have adapted this system to demonstrate the ability of our dyes to undergo bioorthogonal labeling reaction, efficiently. Prior to labeling, Chinese hamster ovary (CHO) cells were incubated with these MegaStokesdyes to test the possible cytostatic effects of these dyes. Experimental results have shown that neither dye is toxic up to 50 µM concentration.


Properties :

- N-Hydroxysuccinimide (NHS) activated acid
- Large Stokes-shift
- excellent photostability
- positively charged
- exceptional in FRET applications
- NIR emission



Absorption / emission maxima :
Molar absorbance :

480 / 600 nm (in ethanol)
47 000 M⁻¹cm⁻¹

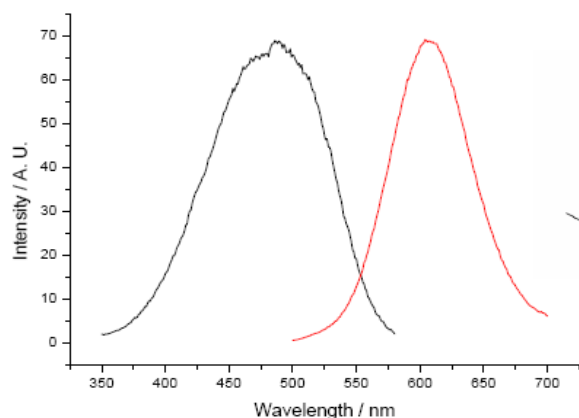

Storage :
Store at -20°C

Precautions and Disclaimer:

This product is for R&D use only, not for drug, household, or other uses. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.

The vibrant M and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. Detailed information on trademarks is available via publicly accessible resources.
© 2018 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

The life science business of Merck KGaA, Darmstadt, Germany operates as MilliporeSigma in the US and Canada.

