




SIGMA QUALITY CONTROL TEST PROCEDURE

Product Information

Enzymatic Assay of L-ALANINE DEHYDROGENASE (EC 1.4.1.1)

PRINCIPLE:

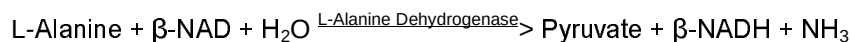

Abbreviations used:

β -NAD = β -Nicotinamide Adenine Dinucleotide, Oxidized Form

β -NADH = β -Nicotinamide Adenine Dinucleotide, Reduced Form

CONDITIONS: T = 25°C, pH = 10.0, $A_{340\text{nm}}$, Light Path = 1 cm

METHOD: Continuous Spectrophotometric Rate Determination

REAGENTS:

- A. 50 mM Sodium Bicarbonate Solution (NaHCO_3)
(Prepare 100 ml in deionized water using Sodium Bicarbonate, Prod. No. S-8875.)
- B. 50 mM Sodium Carbonate Solution (Na_2CO_3)
(Prepare 150 ml in deionized water using Sodium Carbonate, Anhydrous, Prod. No. S-2127.)
- C. 50 mM Sodium Carbonate Buffer, pH 10.0 at 25°C
(Prepare 200 ml by adding 100 ml of Reagent B to 100 ml of Reagent A. Adjust to pH 10.0 at 25°C with Reagent B.)
- D. 500 mM L-Alanine Solution (L-ALA)
(Prepare 1.0 ml in deionized water using L-Alanine, Prod. No. A-7627.)
- E. 30 mM β -Nicotinamide Adenine Dinucleotide, Oxidized Form Solution (β -NAD)
(Dissolve the contents of one 20 mg vial of β -Nicotinamide Adenine Dinucleotide, Stock No. 260-120, in the appropriate volume of deionized water. **PREPARE FRESH.**)
- F. L-Alanine Dehydrogenase Enzyme Solution
(Immediately before use, prepare a solution containing 0.2 - 0.5 unit/ml of L-Alanine Dehydrogenase in cold Reagent C.)

Enzymatic Assay of L-ALANINE DEHYDROGENASE (EC 1.4.1.1)

PROCEDURE:

Pipette (in milliliters) the following reagents into suitable cuvettes:

	<u>Test</u>	<u>Blank</u>
Reagent C (Buffer)	2.70	2.70
Reagent D (L-ALA)	0.10	0.10
Reagent E (β -NAD) 0.10	0.10	

Mix by inversion and equilibrate to 25°C. Monitor the $A_{340\text{nm}}$ until constant, using a suitably thermostatted spectrophotometer. Then add:

Reagent C (Buffer)	-----	0.10
Reagent F (Enzyme Solution)	0.10	-----

Immediately mix by inversion and record the increase in $A_{340\text{nm}}$ for approximately 5 minutes. Obtain the $\Delta A_{340\text{nm}}/\text{minute}$ using the maximum linear rate for both the Test and Blank.

CALCULATIONS:

$$\text{Units/ml enzyme} = \frac{(\Delta A_{340\text{nm}}/\text{min Test} - \Delta A_{340\text{nm}}/\text{min Blank})(3)(\text{df})}{(6.22)(0.1)}$$

3 = Volume (in milliliters) of assay

df = Dilution factor

6.22 = Millimolar extinction coefficient of β -NADH at 340 nm

0.1 = Volume (in milliliter) of enzyme used

$$\text{Units/mg solid} = \frac{\text{units/ml enzyme}}{\text{mg solid/ml enzyme}}$$

$$\text{Units/mg protein} = \frac{\text{units/ml enzyme}}{\text{mg protein/ml enzyme}}$$

UNIT DEFINITION:

One unit will convert 1.0 μmole of L-alanine to pyruvate and NH_3 per minute at pH 10.0 at 25°C.

**Enzymatic Assay of L-ALANINE DEHYDROGENASE
(EC 1.4.1.1)**

FINAL ASSAY CONCENTRATION:

In a 3.00 ml reaction mix, the final concentrations are 47 mM sodium bicarbonate, 17 mM L-alanine, 1.0 mM β -nicotinamide adenine dinucleotide and 0.02 - 0.05 unit of L-alanine dehydrogenase.

REFERENCE:

Bergmeyer, H.U. (1983) *Methods of Enzymatic Analysis*, 2nd edition, Volume I, 427.

NOTES:

1. This assay is based on the cited reference.
2. Where Sigma Product or Stock numbers are specified, equivalent reagents may be substituted.

Sigma warrants that the above procedure information is currently utilized at Sigma and that all Sigma-Aldrich, Inc. products conform to the information in this and other Sigma-Aldrich, Inc. publications. Purchaser must determine the suitability of the information and product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packaging slip.