3050 Spruce Street, St. Louis, MO 63103 USA
Tel: (800) 521-8956 (314) 771-5765 Fax: (800) 325-5052 (314) 771-5757
email: techservice@sial.com sigma-aldrich.com

Product Information

BMPR2 (174-end), active, His-tagged, human PRECISIO® Kinase recombinant, expressed in *Sf*9 cells

Catalog Number **SRP5209** Storage Temperature –70 °C

Synonyms: BMPR-II, BMPR3, BMR2, BRK-3, FLJ41585, FLJ76945, PPH1, T-ALK

Product Description

BMPR2 is a member of the bone morphogenetic protein (BMP) receptor family of transmembrane serine/threonine kinases. The ligands of this receptor are BMPs that are involved in endochondral bone formation and embryogenesis. The loss of interaction and lack of phosphorylation of TCTEL1 by BMPR2 may contribute to the pathogenesis of primary pulmonary hypertension (PPH) and BMPR2 also plays an essential role in human T-cell differentiation.

Recombinant human BMPR2 (174-end) was expressed using baculovirus in *Sf9* insect cells using an N-terminal His tag. The gene accession number is NM_001204. Recombinant protein stored in 50 mM sodium phosphate, pH 7.0, 300 mM NaCl, 150 mM imidazole, 0.1 mM PMSF, 0.25 mM DTT, and 25% glycerol

Molecular mass: ~115 kDa

Purity: 70–95% (SDS-PAGE, see Figure 1)

Specific Activity: 11.9–16.1 nmole/min/mg (see Figure 2)

Precautions and Disclaimer

This product is for R&D use only, not for drug, household, or other uses. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.

Storage/Stability

The product ships on dry ice and storage at -70 °C is recommended. After opening, aliquot into smaller quantities and store at -70 °C. Avoid repeated handling and multiple freeze/thaw cycles.

Figure 1.
SDS-PAGE Gel of Typical Lot 70–95% (densitometry)


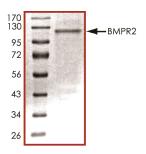


Figure 2.Specific Activity of Typical Lot 11.9–16.1 nmole/min/mg

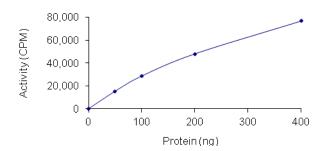

Procedure

Preparation Instructions

Kinase Assay Buffer – 25 mM MOPS, pH 7. 2, 12.5 mM glycerol 2-phosphate, 25 mM MgC1₂, 5 mM EGTA, 2 mM EDTA. Just prior to use, add DTT to a final concentration of 0.25 mM.

Kinase Dilution Buffer – Dilute the Kinase Assay Buffer 5-fold with a 50 ng/µl BSA solution.

Kinase Solution – Dilute the active BMPR2 ($0.1 \,\mu g/\mu l$) with Kinase Dilution Buffer to the desired concentration. Note: The specific activity plot may be used as a guideline (see Figure 2). It is recommended the researcher perform a serial dilution of active BMPR2 kinase for optimal results.

10 mM ATP Stock Solution – Dissolve 55 mg of ATP in 10 ml of Kinase Assay Buffer. Store in 200 μ l aliquots at –20 °C.

 γ -³³P-ATP Assay Cocktail (250 μM) – Combine 5.75 ml of Kinase Assay Buffer, 150 μl of 10 mM ATP Stock Solution, 100 μl of γ -³³P-ATP (1 mCi/100 μl). Store in 1 ml aliquots at –20 °C.

Substrate Solution – Myelin Basic Protein (MBP) substrate diluted in distilled water to a final concentration of 1 mg/ml.

1% phosphoric acid solution – Dilute 10 ml of concentrated phosphoric acid to a final volume of 1 L with water.

Kinase Assay

This assay involves the use of the ³³P radioisotope. All institutional guidelines regarding the use of radioisotopes should be followed.

- 1. Thaw the active BMPR2, Kinase Assay Buffer, Substrate Solution, and Kinase Dilution Buffer on ice. The γ -³³P-ATP Assay Cocktail may be thawed at room temperature.
- In a pre-cooled microcentrifuge tube, add the following solutions to a volume of 20 μl:

10 μ l of Kinase Solution

5 μl of Substrate Solution

5 μl of cold water (4 °C)

- 3. Set up a blank control as outlined in step 2, substituting 5 μ l of cold water (4 °C) for the Substrate Solution.
- 4. Initiate each reaction with the addition of 5 μ l of the γ - 33 P-ATP Assay Cocktail, bringing the final reaction volume to 25 μ l. Incubate the mixture in a water bath at 30 °C for 15 minutes.
- After the 15 minute incubation, stop the reaction by spotting 20 μl of the reaction mixture onto an individually precut strip of phosphocellulose P81 paper.

- Air dry the precut P81 strip and sequentially wash in the 1% phosphoric acid solution with constant gentle stirring. It is recommended the strips be washed a total of 3 times of ~10 minutes each.
- 7. Set up a radioactive control to measure the total γ - 33 P-ATP counts introduced into the reaction. Spot 5 μ l of the γ - 33 P-ATP Assay Cocktail on a precut P81 strip. Dry the sample for 2 minutes and read the counts. Do not wash this sample.
- 8. Count the radioactivity on the P81 paper in the presence of scintillation fluid in a scintillation counter.
- 9. Determine the corrected cpm by subtracting the blank control value (see step 3) from each sample and calculate the kinase specific activity

Calculations:

1. Specific Radioactivity (SR) of ATP (cpm/nmole)

SR = $\frac{\text{cpm of 5} \ \mu \text{l of } \ \gamma^{-33} \text{P-ATP Assay Cocktail}}{\text{nmole of ATP}}$ cpm – value from control (step 7)
nmole – 1.25 nmole (5 μ l of 250 μ M ATP Assay Cocktail)

2. Specific Kinase Activity (SA) (nmole/min/mg)

nmole/min/mg =
$$\Delta$$
cpm × (25/20)
SR × E × T

SR = specific radioactivity of the ATP (cpm/nmole ATP) Δ cpm = cpm of the sample – cpm of the blank (step 3) 25 = total reaction volume

20 = spot volume

T = reaction time (minutes)

E = amount of enzyme (mg)

References

- Machado, R.D. et al., Functional interaction between BMPR-II and Tctex-1, a light chain of dynein, is isoform-specific and disrupted by mutations underlying primary pulmonary hypertension. Hum. Molec. Genet., 12, 3277-3286 (2003).
- Cejalvo, T. et al., Bone morphogenetic protein-2/4 signalling pathway components are expressed in the human thymus and inhibit early T-cell development. Immunology, 121, 94-104 (2007).

PRECISIO is a registered trademark of Sigma-Aldrich Co. LLC.

RC,MAM 10/11-1