

User Guide

DirectLoad™ Plus 1 Kb DNA Ladder

DPLUS-1K**500 µL/vial**

Storage

-20 °C for up to 24 months, 4 °C for up to 12 months, or room temperature for up to 6 months.

Description


PCR products and double-stranded DNA digested with appropriate restriction enzymes equilibrated to 10 mM Tris-HCl (pH 8.0) and 10mM EDTA after phenol extraction to yield 12 fragments, suitable for use as molecular weight standards for agarose gel electrophoresis. The DNA includes fragments ranging from 100-10,000 base pairs. The 1,000 and 3,000 base pair bands have increased intensity to serve as reference points.

Note: The approximate mass of DNA in each band is supplied for estimating the mass of DNA in comparably intense samples of similar size.

Range	100-10,000 bp
Number of bands	13
Concentration	100 µg/mL
Tracking dye	Bromophenol Blue and Xylene Cyanol FF

This product is only suitable for rough estimation of nucleic acid molecular weight.

Suggested Load: 5 µL/well

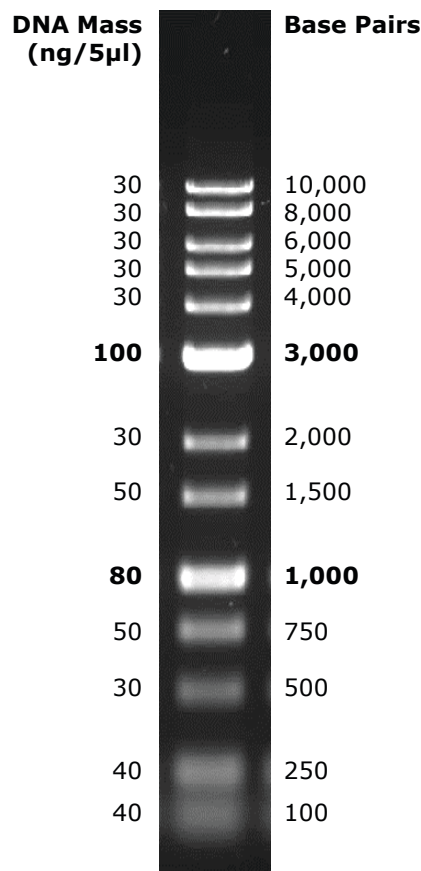

1% TAE agarose gel

Product Ordering

Purchase online at SigmaAldrich.com.

Description	Qty	Catalogue Number
DirectLoad™ Plus 50 bp DNA Ladder	500 µL vial	DPLUS-50
DirectLoad™ Plus 100 bp DNA Ladder	500 µL vial	DPLUS-100
DirectLoad™ Plus 100 bp High DNA Ladder	500 µL vial	DPLUS-100H
DirectLoad™ Plus 100 bp High DNA Ladder with fluorescent DNA stains	600 µL vial	DPLUS-100HS
DirectLoad™ Plus 1 Kb DNA Ladder	500 µL vial	DPLUS-1K
DirectLoad™ Plus 1 Kb DNA Ladder with fluorescent DNA stains	600 µL vial	DPLUS-1KS
BioReagent Ethidium Bromide Solution	10 mL	E1510-10ML
BioReagent Agarose, low electroendosmosis (EEO)	50 g	A9539-50G
pHast Pack™ Tris-Acetate EDTA Buffer, pH 8.3	20 pouches	PPB008-20PAK
pHast Pack™ Tris-Borate EDTA Buffer, pH 8.3	20 pouches	PPB009-20PAK
pHast Pack™ Tris-Borate EDTA w/ 1% agarose	20 pouches	PPB012-20PAK
pHast Pack™ Tris-Acetate EDTA w/ 1% agarose	20 pouches	PPB013-20PAK
GelRed® Nucleic Acid Stain 3X Water	4 L	SCT121
GelRed® Nucleic Acid Stain 10000X DMSO	0.5 mL	SCT122
GelRed® Nucleic Acid Stain 10000X Water	0.5 mL	SCT123
GelGreen® Nucleic Acid Stain 10000X DMSO	0.5 mL	SCT124
GelGreen® Nucleic Acid Stain 10000X Water	0.5 mL	SCT125
SYBR® Green I Nucleic Acid Gel Stain	0.5 mL 1 mL	S9430-.5ML S9430-1ML

Precautions and Disclaimer

This product is for research use only. Please consult the Safety Data Sheet for information regarding hazards and safe handling practices.

Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The information in this document is subject to change without notice and should not be construed as a commitment by the manufacturing or selling entity, or an affiliate. We assume no responsibility for any errors that may appear in this document.

Contact Information

For the location of the office nearest you, go to SigmaAldrich.com/offices.

Technical Assistance

Visit the tech service page on our web site at SigmaAldrich.com/techservice.

Standard Warranty

The applicable warranty for the products listed in this publication may be found at SigmaAldrich.com/terms.

Merck, DirectLoad, pHast Pack and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.
© 2022 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

The life science business of Merck operates as MilliporeSigma in the U.S. and Canada.

20667263w Rev 03/22

