

Technical Data Sheet

CE SABOURAUD-2 % Dextrose Agar

Ordering number: 1.07315.0500

The medium was recommended by JANKE (1961) for the cultivation of dermatophytes. GEORG et al. (1954) recommended the addition of cycloheximide, penicillin and streptomycin for the inhibition of non-pathogenic accompanying bacteria.

This medium is used for the cultivation of fungi (yeasts, molds), especially recommended for the fungi associated with skin infections. SABOURAUD-2 % Dextrose Agar is often used with antibiotics for the isolation of pathogenic fungi from material containing a variety of other fungi or bacteria. It is only slightly selective against bacteria.

IVD in vitro diagnosticum - For professional use only

Typical Composition

Peptone	10.0 g/l
D(+)-glucose	20.0 g/l
agar-agar	17.0 g/l

Preparation

Suspend 47 g/litre, autoclave (15 min at 121 °C).

Do not overheat.

pH: 5.6 ± 0.2 at 25 °C.

The plates are clear and yellowish-brown.

Experimental Procedure and Evaluation

The plates are inoculated with sample material according to the instructions. The fungi colonies which have grown are judged macro- and microscopically.

Incubation for up to 7 days at 28 °C aerobically.

For preparation of version acc. to EMMONS adjust pH to 6.9 ± 0.2 . Chloramphenicol (50 mg/l) may be added in addition.

Storage

Usable up to the expiry date when stored dry and tightly closed at +15 to +25° C. Protect from light.

After first opening of the bottle the content can be used up to the expiry date when stored dry and tightly closed at +15 to +25 °C.

Specimen


e.g. Skin, nails, sputum, exdates, open lesions.

Clinical specimen collection, handling and processing, see general instructions of use.

Quality Control

Control Strains	ATCC #	Incubation	Expected Results
<i>Trichophyton mentagrophytes</i>	18748	7 d at 28 °C	Growth good / very good
<i>Trichophyton rubrum</i>	28188	7 d at 28 °C	Growth fair / good
<i>Microsporum gallinae</i>	12108	7 d at 28 °C	Growth good / very good
<i>Trichophyton ajelloi</i>	28454	7 d at 28 °C	Growth fair / very good
<i>Microsporum canis</i>	36299	7 d at 28 °C	Growth good / very good
<i>Geotrichum candidum</i>	DSMZ 1240	7 d at 28 °C	Growth good / very good
<i>Candida albicans</i>	10231	7 d at 28 °C	Growth good / very good
<i>Aspergillus niger</i>	16404	7 d at 28 °C	Growth good / very good
<i>Penicillium commune</i>	10428	7 d at 28 °C	Growth good / very good

Please refer to the actual batch related Certificate of Analysis.

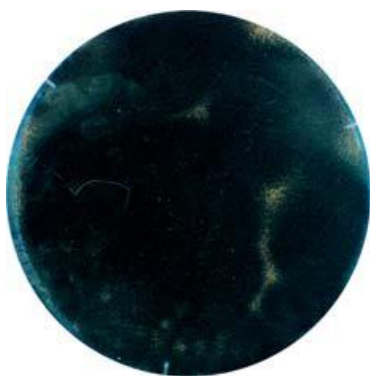

Aspergillus niger ATCC 16404



Geotrichum candidum DSMZ 1240

Literature

GEORG, L.K., AJELLO, L., a. PAPAGEORGE, C.: Use of cycloheximide in the selective isolation of fungi pathogenic to man. - J. Lab. Clin. Med., 44; 422-428 (1954).

JANKE, D.: Pilznährboden nach SABOURAUD, modifiziert MERCK, ein neuer Trockennährboden zur Züchtung von Dermatophyten. - Zschr. Haut- u. Geschl.-Krankh., 15; 188-193 (1961).



We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

Merck Millipore, the M logo and Merckoplate® are registered trademarks of Merck KGaA, Darmstadt, Germany. ATCC is a registered trademark of ATCC, Manassas, VA, USA. Lit. No. TN1468EN00

Ordering Information

Product	Cat. No.	Pack size
SABOURAUD-2 % Dextrose Agar	1.07315.0500	500 g

Merck KGaA, 64271 Darmstadt, Germany
Fax: +49 (0) 61 51 / 72-60 80
mibio@merckgroup.com
www.merckmillipore.com/biomonitoring

Find contact information for your country at:
www.merckmillipore.com/offices
For Technical Service, please visit:
www.merckmillipore.com/techservice



We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

Merck Millipore, the M logo and Merckoplate® are registered trademarks of Merck KGaA, Darmstadt, Germany. ATCC is a registered trademark of ATCC, Manassas, VA, USA. Lit. No. TN1468EN00